User Manual

Worklist Server User Manual_en

		Page 1 of 9
Version: 1.2	Reviewed by: Veronica Simonet	Approved by: Rafael Sanguinetti
	Date: 28/11/2019	Date: 28/11/2019

User Manual
Worklist Server

User Manual

Worklist Server Manual_en

Version 1.2 Page 2 of 9

Table of Content

1.	Intr	roduction	3			
2.	Computer Platform					
		ftware installation				
	WLS Functionality					
		Main Window				
		Setup				
		2.1 Server Configuration				
		2.2 Delete Rules				
4		WEB Interface				

User Manual

Worklist Server Manual en

Version 1.2 Page 3 of 9

1. Introduction

This application is a standalone **Worklist Server**, integrating scheduled studies with medical imaging equipment. A Worklist Server eliminates manual entry errors, speeding up study completion.

Scheduled studies can be sent by HL7 ORM commands or using HTTP protocol Post/Get calls. Multiple modalities can be connected, guerying simultaneously the worklist server.

The list can be updated continuously, deleting completed studies, using a MPPS simulation mode.

2. Computer Platform

- Windows 7, Windows 8 or Windows 10, (32 or 64 bits)
- Intel i3, i5 or i7 processor.
- 4 GBytes RAM.
- 300 GBytes HDD.

3. Software installation

Software installation:

- 1. Download the Installer from http://charruasoft.com/products/wls/index.php
- 2. Unzip the ZIP file downloaded.
- 3. Run WLS install.exe
- 4. Follow the Wizard to correctly install the software.

Note: It is suggested to install the Software in the path recommended in the installation Wizard (C:\CharruaSoft).

4. Worklist Server Functionality

4.1 Main Window



Figure 1. Main Window

After the installation is completed, the main window will appear. You must go to the Setup tab to configure the Worklist Server.

4.2 Setup

4.2.1 Server Configuration

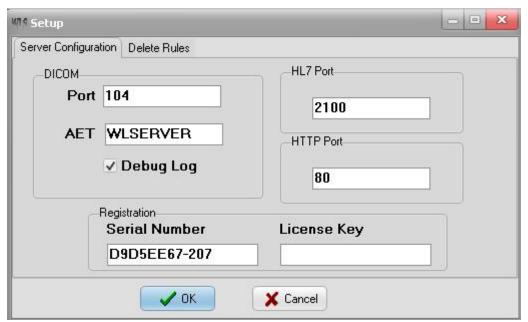


Figure 2. Server Configuration Window

The HL7 Listener uses HL7 Port.

The HTTP Port is for the WEB server. The WEB server is for configuration and for adding or querying the Worklist database.

The DICOM Worklist Server configuration needs the AET and Port number.

4.2.2 Discontinued and Completed Rules

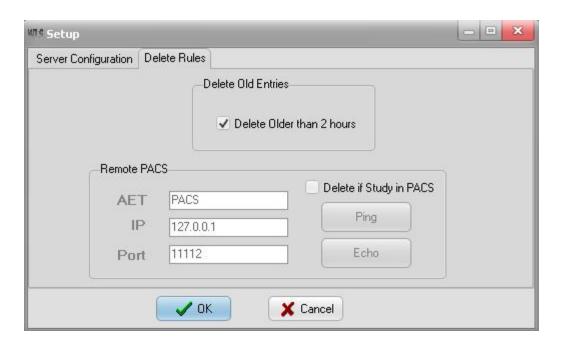


Figure 2. Delete Rules

Remote PACS configuration information (AET, IP Address and Port number), are needed in case of using the MPPS simulation. This will delete worklist entries once the study images arrive at the PACS system.

The other option for automatically deleting the entries is selecting the "Delete Older than 2 hours" checkbox.

User Manual

Worklist Server Manual_en

Version 1.2 Page 7 of 9

4.3 WEB Interface

	Worklist WEB Server - admin	
	Manage Users View and manage User database	
■	Query Worklist Browse worklist entries	
■	Add to Worklist Add a new worklist entry	


This is the administrative view, it features:

- User management; add, query, modify or delete users.
- Query Worklist: query, modify or delete Worklist Entries.
- Add to Worklist: add new entries into the Worklist.


User Manual

Worklist Server Manual_en

Version 1.2 Page 8 of 9

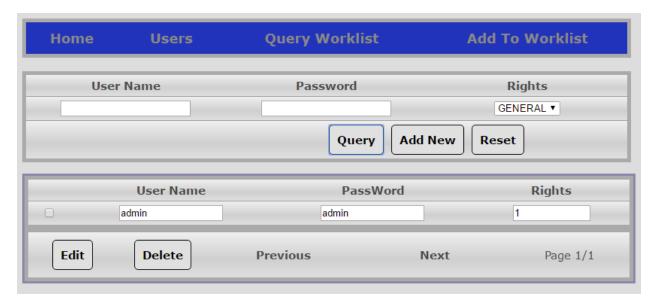

User Management

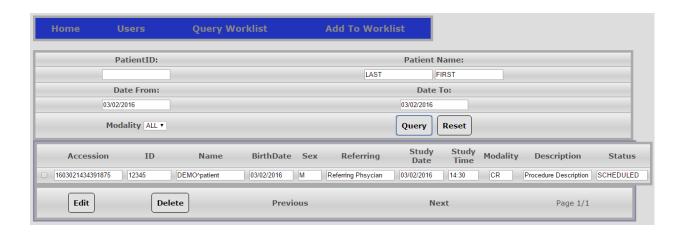


• Query Worklist; filter the worklist and list scheduled studies.

User Manual

Worklist Server Manual_en

Version 1.2 Page 9 of 9

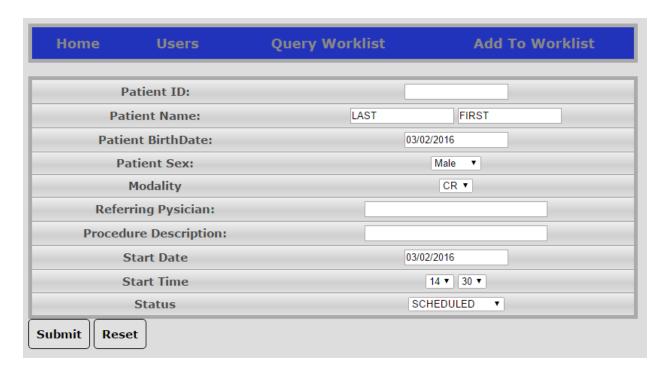

Add to Worklist, add a new study to the list.

General users only can query or add to the Worklist.

8. Document Change History

Document Version	Date	Changed performed by	Change description
1.0	1/3/2017	Rafael Sanguinetti	Document Creation
1.1	02/03/2016	Rafael Sanguinetti	Small updates were made.
1.2	28/11/2019	Rafael Sanguinetti	The chapters referring physicians and procedure descriprion data base were removed